



INVACARE® TOP END®  
**SCHULTE 7000 SERIES**  
BASKETBALL WHEELCHAIR




The **Invacare® Top End® Schulte 7000 Series Basketball Wheelchair** is faster, lighter and stronger. In 2011 it hit the market with new material, enhanced headtube design and axle securing system, reducing overall weight by two pounds and improving performance. Each Schulte 7000 Series Wheelchair is custom built to each individual athlete's specifications and is overseen by Paul Schulte who is the design engineer for Invacare Top End and three-time Paralympian as a member of the USA National Men's Wheelchair Basketball Team.

## Features

Custom-built fixed frame in 7005 aluminum  
Adjustable height fixed swivel anti-tip 58mm fifth wheel  
Fixed back with fixed or adjustable height  
Adjustable tension seat and back upholstery  
Fixed axle position with fixed wheelbase  
Fixed camber: 15°, 18° or 20°  
High performance Sun® wheels: 24", 25", 26" or 700C  
Four bed offensive wing  
Anodized pushrims

High pressure clincher tires  
Precision fork caster system  
72mm casters  
Adjustable height tubular footrest  
Powder coat paint  
Rear frame protector  
Axle securing system  
250 lb. weight capacity

## Model #: PS7

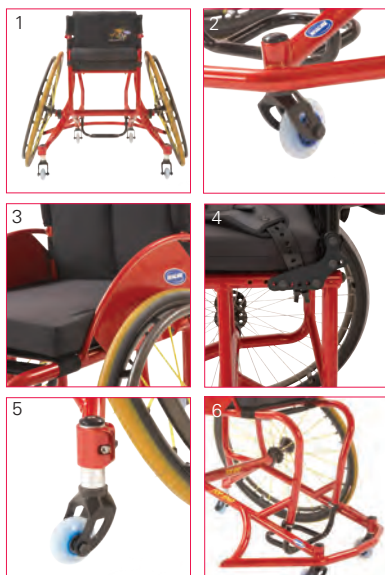

1. Double swivel anti-tip/ fifth wheel option
2. Precision fork caster system
3. Optional welded fixed tube sideguard
4. Adjustable center-of-gravity option
5. Adjustable height fixed swivel anti-tip
6. 7005 Aluminum heat treating process



**Invacare Corporation**  
[www.invacare.com](http://www.invacare.com)

### USA

One Invacare Way  
Elyria, Ohio  
44035-4190  
(800) 333-6900

### Canada

570 Matheson Blvd. E.,  
Unit 8  
Mississauga, Ontario  
L4Z 4G4 Canada  
(800) 668-5324

### Top End

[www.topendwheelchair.com](http://www.topendwheelchair.com)  
Sports and Recreation Products  
4501 63rd Circle North  
Pinellas Park, Florida  
33781-5914  
(727) 522-8677  
(800) 532-8677  
Fax: (727) 522-1007

©2015

Invacare Corporation.  
All rights reserved.  
Trademarks are identified  
by the symbols ™ and ®.  
All trademarks are owned  
by or licensed to Invacare  
Corporation unless  
otherwise noted. Spinerigy  
is a registered trademark of  
Spinerigy Inc. Sun and Equalizer  
are registered trademarks  
of Sun Components.  
Specifications are subject to  
change without notification.  
Form No. 12-480  
Rev. 03/16  
151172

## Options

Adjustable Center-of-Gravity  
Spinerigy® SLX wheels  
Sun® Equalizer™ wheels  
12 Tab pushrims  
Tapered seat width  
Welded seat rigidizer bars with solid seat  
Double swivel anti-tip/fifth wheel  
Click straps  
Double amputee front-end

Recessed platform footrest  
Adjustable angle footrest  
Clamp-on fixed aluminum side guards  
Fixed side guards with welded tube  
Fixed aluminum or carbon fiber side guards  
Titanium axles  
Titanium pushrims  
Front frame protectors  
Alignment gauge

