



Domus 4


Alternating Pressure Redistribution System

- Suitable for patients with medium to high risk of pressure injury
- Simple-to-use weight setting function
- CPR knob for quick deflation & maximum security
- Heel relief function for optimal heel protection



Care for a Healthy Life

.PEX

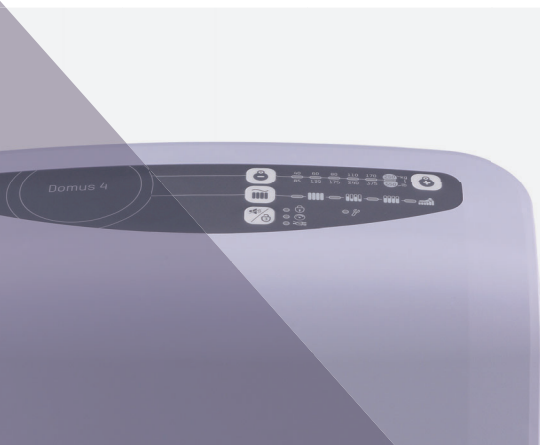

Domus 4

Alternating Pressure Redistribution System



Cell-on-Cell Design
The bottom cells prevent the patient from bottoming out in case of power failure



Micro Low Air Loss
Helps wick away moisture from the surface to keep patient's skin less susceptible to injuries



Heel Relief Function
Individual detachable air cells prevent pressure injuries on the heels



Visual and Audible Alarm

- Low Pressure
- Power Failure
- Tech Support



CPR Knob
Designed for quick mattress deflation in case of emergency and prevents accidental triggering



4-Way Stretch Cover
Low friction and shear cover with high moisture-vapour transmission rate for better microclimate management



Biocompatibility



Vapour Permeable



Water Resistance



Flame Resistance




Enhanced Cross-Infection Control
Made with all-around-zipper and welded seams for easy removal and prevents fluid ingress



Quick Connector with Cap
Facilitates system installation and seals the mattress for patient transportation up to 8 hours

Pressure Relief During Alternating

The Domus 4 is designed to maximize patient's safety and offers multiple therapy modes that are capable of fulfilling various clinical needs for better patient care.

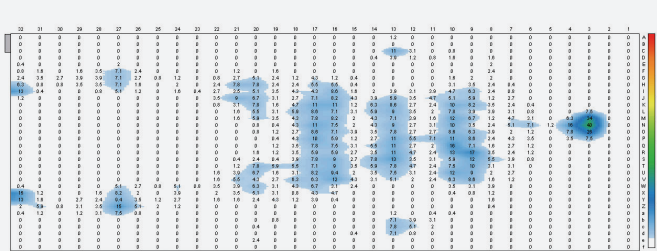

5 Modes to Fulfill Individual Needs

- Max Firm
- Alternating
- Continuous Low Pressure (CLP)
- Seat Inflation w/ Alternating
- Seat Inflation w/ CLP



Alternating

- Alternately inflates and deflates air cells allowing pressure redistribution and encourage tissue reperfusion



Max Firm

- A static firm surface for nursing procedures. It will automatically return to alternating mode after 20 minutes



Continuous Low Pressure (CLP)

- Provides a constant static surface in lower pressure with a greater contact area for patients to be immersed and enveloped better



Seat Inflation

- Adjust the pressure to suit the patient in a sitting position. It can be used with alternating or CLP mode

Specifications

Model	Domus 4
Pump	Dimension 29 x18.5 x 12.6 cm
	Weight 2.2 kg
	Case Material Flame Retardant ABS
	Supply Voltage AC 220 – 240, 50 Hz
	Operating Cycle 10 mins
Mattress	8" Replacement
	Dimension 200 x 80/85/90 x 20.3 cm
	Cell Height 8" x 20 Cells
	Weight 6.9 kg
	Cover Material 4-Way Stretch
	Cell Material Nylon TPU
	Maximum Patient Weight 250 kg
	Flame Retardant Standards EN 591-1; EN 591-2

Specifications are subject to change for improvement without notice

